

# Shop the Way You Want



Feel Before You Buy

Komet's reps carry bur wallets to every office to show you the product before ordering.



Shop on Your Schedule

Komet's website offers Komet's entire portfolio of premium products at your convenience.

# **High-Performance Ceramics**

#### **Shop Now**

# **Trimming DCBs**

- Effectively, gently and intuitively trim high-performance ceramics
- Only light pressure needed due to interspersed diamonds and sintered ceramic bond



Extended shapes and sizes now available!

# Adjusting **ZR** Diamonds

- Reduce material rapidly with little effort due to extremely durable diamonds
- Feature densely packed diamond layer with permanently bonded, high-quality diamond particles





### Diamond Discs

#### 6942

- Sintered disc used for separating ceramics
- Long service life due to diamond interspersed edge

### 6924

- Spiral reinforced disc
- Increases visibility when separating and contouring ceramics






# Polishing

6942

### 7R Flash Polishers™

High performance, two-step polishers interspersed with diamonds.


- Spiral polishers: For occlusal surfaces
- Wheel polishers: For lingual or buccal mass polishing
- Flame polishers: For buccal polishing

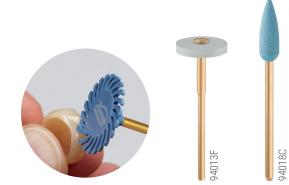



# Extra-Oral All-Ceramic Adjustment Kit

Comprehensive kit with coarse polishers to efficiently adjust and polish high-performance ceramics.

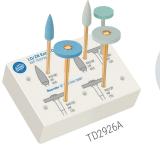



Includes a diamond ceramic abrasive that proficiently cuts while generating minimal heat.

DCB4

# **Acrylics/Resin**

# Separating

### 987P

- Large, serrated diamond disc used for separating resin models

# Can be used in a clockwise or counterclockwise rotation Adjusting Carbide Cutters

# Coarse Toothing

- ACR toothing resists jamming and is easily controlled due to its staggered toothing, making it ideal for denture acrylics
- GSQ toothing provides easy and controlled surface shaping, effective cutting and optimal surface quality for soft and temporary acrylics

### Medium Toothing

E toothing removes excessive material and trims acrylic gently with little pressure

# Fine Toothing

EF toothing provides super precise finishing and short granular chips ensure smooth cutting









H79GSQ



H79ACR



# Polishing

### 9486

Abrasive buff embedded in silicone to prevent shedding while removing, rounding and smoothing plastics.

### 9489C/M


Coarse and fine grit polishers with improved durability for pre-polishing.

# 9644F

Used for high-shine polishing.

# Brush Mop Polishers

- Coarse, medium and fine grits for polishing acrylics
- Used for smoothing, buffing or polishing specific spots on clear aligners that feel uncomfortable to patients

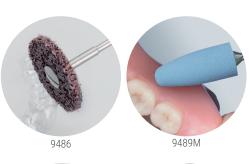





# **Stone and Model**

# Cutting/Separating



**Shop Now** 



987P

- Serrated diamond disc used for cutting plaster models
- Can be used in clockwise or counterclockwise rotation



## Carbide Cutters

### Coarse Toothing

- · SGFA toothing reduces bulk material on wet plaster
- GEA toothing reduces bulk material on dry plaster

# Medium Toothing

 E toothing removes excessive material and trims dry plaster

Features staggered toothing that divides the instrument blades into individual, offset cutting segments.

## Fine Toothing

 EF toothing provides super-precise finishing and short granular chips ensure smooth cutting



6924

- Spiral reinforced disc that provides efficient material reduction and greater stability
- Features honeycomb design and extra-large chip spaces



H79SGFA



H251GEA



H79E



H73EF